www.roshantebsepanta.com

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Stand by | DUAL | Req./array buffer 2/2 | ALL TESTS | 12:39:15 | 21/01/2010

3*

Request samples

Man menu
Special functions
Rotor/System
Installation
Programming
Service
Funct. ck/adj.
Error history
Quality control
Evaluate samples
Blank/Calib info
Eval. disk samp.
Request samples
Sample handling
Load menu
Unload menu
Reagent info

Sample ID:

Patient name:

Date of birth:

Sex:

Physician:

Repeats:

Sample type:

ALBUMIN
 ALKALINE PHOS.
 ALAT/GPT
 ASAT/GOT
 AMYLASE
 TOTAL BILIRUBIN
 DIRECT BILIRUBI
 CALCIUM
 CHOLESTEROL
 HDL CHOLESTEROL
 LDL CHOLESTEROL
 CK TOTAL
 CREAT. ENZ.
 GAMMA GT
 GLUCOSE

Pos	Sample ID	Name
b		Blank
s		ELUCAL
*		

F1 Worklist | F2 STAT | Shift+F3/Alt+F3 Delete/ Delete All | F4 Repeat Mode | F5 Previous Sample | F6 Next Sample | F7 Profile | F8 New Sample | F9 Sample Handling | F10 Main Menu

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Stand by | DUAL | Req./array buffer 2/2 | ALL TESTS | 12:40:19 | 21/01/2010

3*

Request samples

Man menu
Special functions
Rotor/System
Installation
Programming
Service
Funct. ck/adj.
Error history
Quality control
Evaluate samples
Blank/Calib info
Eval. disk samp.
Request samples
Sample handling
Load menu
Unload menu
Reagent info

Control name:

Batch number:

Expiry date:

Repeats:

Sample type:

Pos	Sample ID	Name
b		Blank
s		ELUCAL
*		Control

F1 Worklist | F2 STAT | Shift+F3/Alt+F3 Delete/ Delete All | F4 Repeat Mode | F5 Previous Sample | F6 Next Sample | F8 New Sample | F9 Sample Handling | F10 Main Menu

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
The screenshot shows the software interface with a circular sample layout on the left and a 'Sample Handling' table on the right. The table is circled in green.

Pos	Sample ID	Name
b		Blank
s		ELICAL
c		ELITROL 1
c	999-001	ELITROL 2 CLOONEY G.

At the bottom of the interface, there are function keys: F1 (Print Loadlist), F5 (Reagent Info), F6 (Quality Control), F7 (Evaluate Samples), F8 (Request Samples), and F10 (Main Menu).

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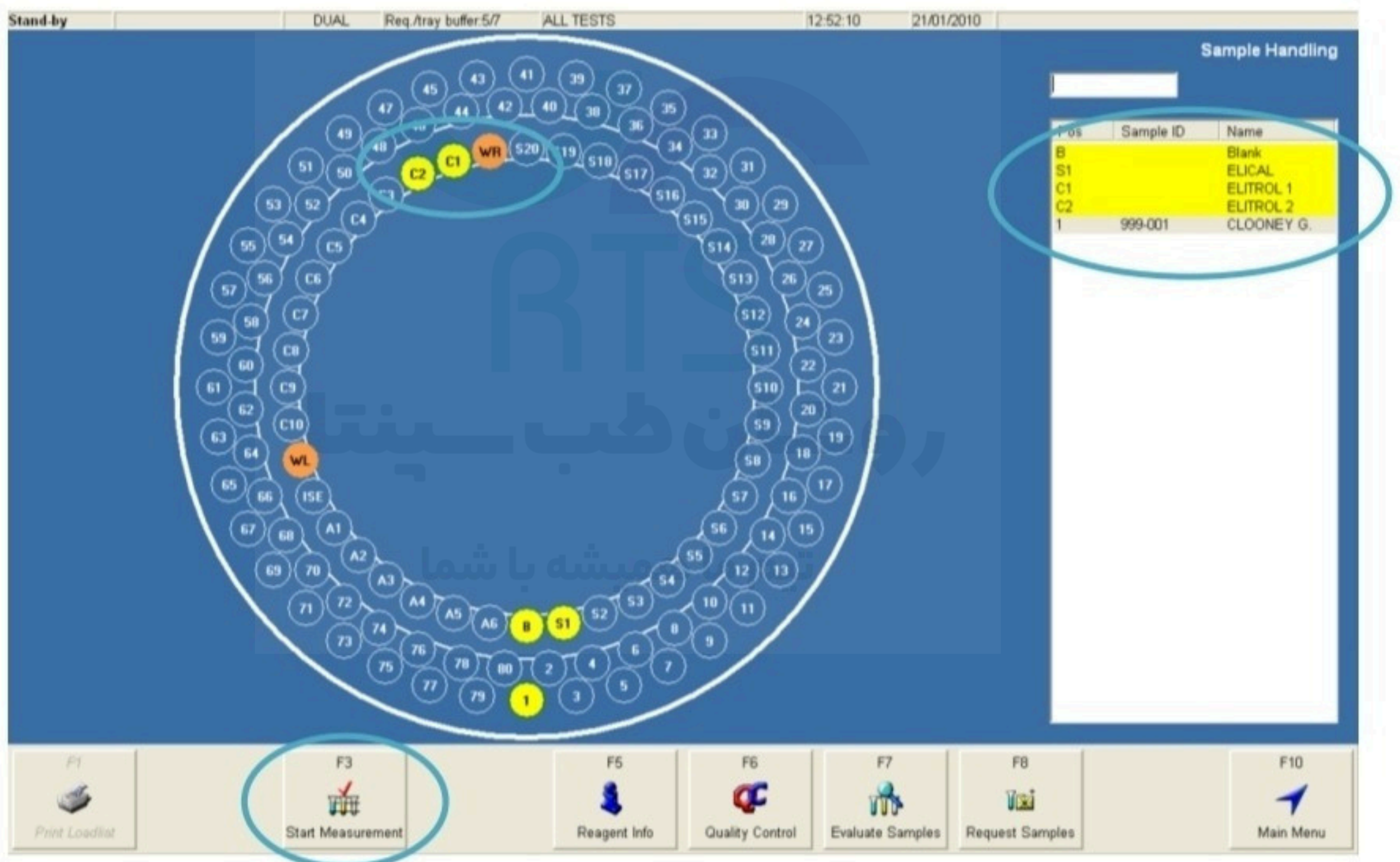Stand by | DUAL | Req /tray buffer 5/7 | ALL TESTS | 12:52:10 | 21/01/2010

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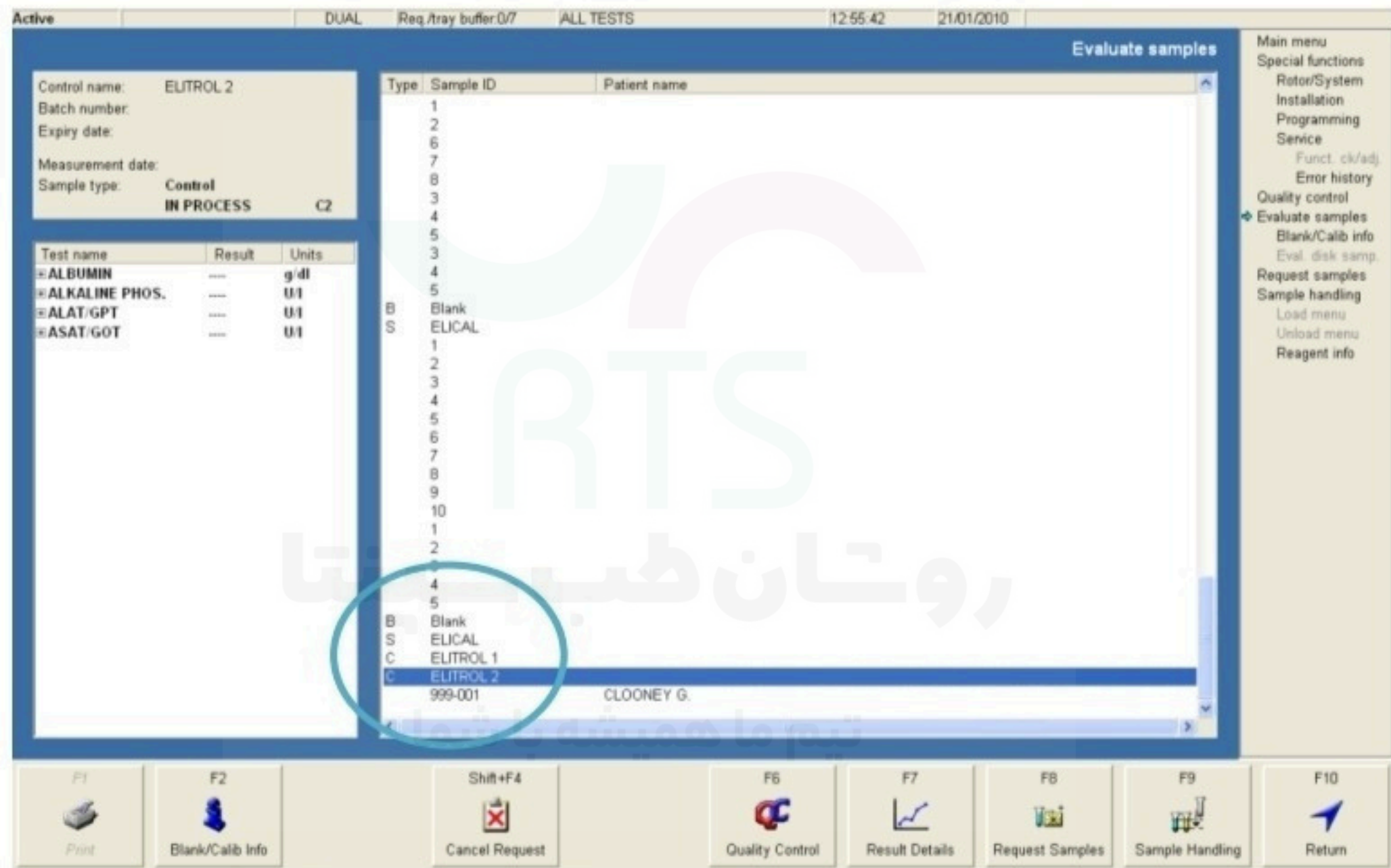
The screenshot displays the 'Evaluate samples' interface. At the top, the status bar shows 'Active', 'DUAL', 'Req./tray buffer:0/7', 'ALL TESTS', '12:55:42', and '21/01/2010'. The left sidebar contains control information for 'EUTROL 2' and a table of test results.

Control name:	EUTROL 2
Batch number:	
Expiry date:	
Measurement date:	
Sample type:	Control IN PROCESS C2

Test name	Result	Units
ALBUMIN	---	g/dl
ALKALINE PHOS.	---	U/l
ALAT/GPT	---	U/l
ASAT/GOT	---	U/l

The central table lists samples with the following data:

Type	Sample ID	Patient name
1		
2		
6		
7		
8		
3		
4		
5		
3		
4		
5		
B	Blank	
S	ELICAL	
1		
2		
3		
4		
5		
6		
7		
8		
9		
10		
1		
2		
3		
4		
5		
B	Blank	
S	ELICAL	
C	ELITROL 1	
C	ELITROL 2	CLOONEY G.
	999-001	

The right sidebar contains a 'Main menu' with options: Special functions, Rotor/System, Installation, Programming, Service, Funct. ck/ed, Error history, Quality control, Evaluate samples (highlighted), Blank/Calib info, Eval. disk samp., Request samples, Sample handling, Load menu, Unload menu, and Reagent info.

The bottom toolbar includes function keys: F1 (Print), F2 (Blank/Calib Info), Shift+F4 (Cancel Request), F6 (Quality Control), F7 (Result Details), F8 (Request Samples), F9 (Sample Handling), and F10 (Return).

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